
User manual

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MITOCHONDRIAK® 1.0

Mitochondriak® Slim

Mitochondriak® Maxi

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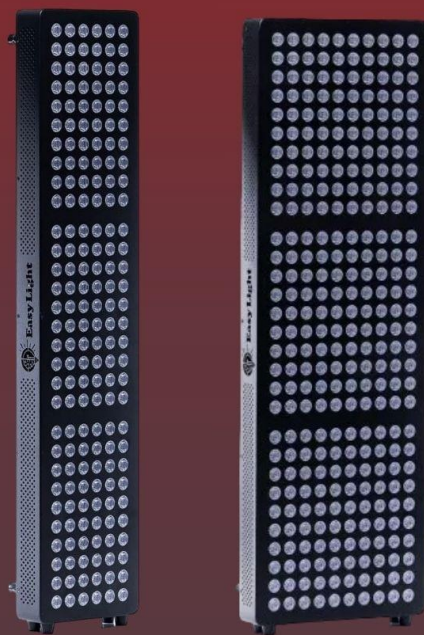@easylight.sk

User manual content:

Easy Light Mitochondriak® panel models	24
Technical parameters	25
Setting the timer	26
Remote control	27
How red light therapy works	28
Installation + hanging of the device	31
Instructions for use	32
FAQ - Frequently asked questions	33
What's in the package?	37
Maintenance and cleaning instructions	37
Warnings and Cautions	38
Warranty	40

Mitochondriak® 1.0 panel models:

Mitochondriak® Maxi

91*30 cm

300 LED pcs

Mitochondriak® Slim

91*20 cm


180 LED pcs

Technical parameters:

Model:	Mitochondriak® Slim	Mitochondriak® Maxi
Dimension:	910x195x60 mm	910x298x65 mm
Intensity:	>120 mW/cm ² from a distance 15 cm	
Flicker:	Reduced (1Hz, 0,9 to 3,5%, Index = 0 až 0,03)	
LED power:	540 W	900 W
Number of LEDs:	180x3W	300x3W
LED illumination angle:	30 degrees	
Wavelengths:	630, 670, 810, 830, 850nm (approx. 20:30:10:10:30)	
Timer control:	5/10/15/20/25/30/ --- 60 minutes	
Input voltage:	AC 110-240 V	
Frequency:	50/60 Hz	
EMF:	at the recommended distance almost zero	
Working temperature:	-20 °C to +40 °C	
Material:	Metal/Plastic	
Weight:	8 Kg	11 Kg
LED lifespan:	50 000 hours (137 years with 60 min./day lighting)	
Warranty:	3 years	
Certificate:	CE, FCC, ROHS,	

Timer setting and control side panel:

5. Turn on the power at the top of the panel, the LED screen on the side display shows “OFF” by default.
6. Set the light you want to turn on. Press the “RED” button if you want only red light to shine from the panel and then the corresponding “RED” light will turn green. If you want the near infrared light to shine, press the “NIR” button. Then the corresponding green NIR light will turn on. If you turn on both, both red and infrared will light up (both RED and NIR lights will be lit).
7. Press the "TIME+" or "TIME-" button and set the timer according to your requirements (5 - 10 - 15 - 20 - 25...30 minutes).
8. The panel will automatically turn itself off after the set time has elapsed and the LED display will show “OFF”.

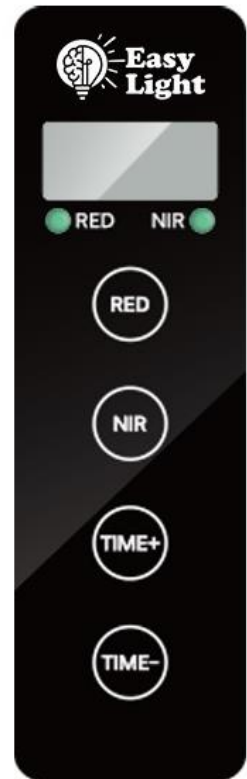Remote control:


For the remote control, use a **23A 12V** battery. You can buy it at a regular electrical store. The battery is not included in the package.

The remote control function is automatically preset and you can use the controller after turning on the panel.

However, sometimes it may happen that the remote control gets unpaired or stops working, or you may want to unpair/pair it again. In this case, please follow these simple steps:

To disable/enable switching using the remote control, simply take the remote control, point it at a device that is plugged into the network, but turned OFF with the main red switch (at the back) to the OFF position. Then press and hold the RED button on the remote control for about 6 to 7 seconds and turn the main switch on the panels to the ON position. That's it.

If this does not work, it is possible that synchronization is disabled on the panel. First, on the side of the switched on panel, press TIME+ for about 6-7 seconds, which will enable the synchronization option.

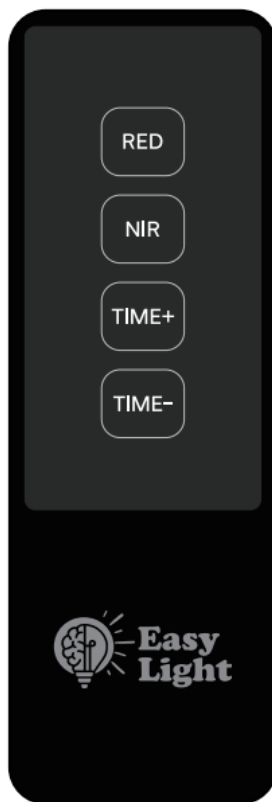Note: The EMF from the remote control is minimal, but as with any wireless electronics, we recommend holding it in your hand only when switching!

How red light therapy works:

Introduction:

Thank you for purchasing an **Easy Light Mitochondriak®** light therapy device. You have chosen panels that have been purposely designed to provide the most optimal experience and effect from red and infrared light therapy.

How it works:


Red and infrared light therapy (also known as RLT – Red Light Therapy, or LLLT – Low Level Laser Therapy) is a therapy that uses red or near-infrared light (so-called NIR).

The LEDs in the panels are designed to emit light only at very specific, pre-programmed, therapeutic wavelengths of light. EasyLight panels utilize wavelengths that have been shown to have positive biological results.

The vast majority of clinical literature suggests that light in the mid-600nm to 850nm range is extremely effective, mainly due to its effect on cellular respiration in mitochondria, but also on water, skin, and collagen. This is why EasyLight Mitochondriak® has selected up to 5 wavelengths that have multiple biological effects.

EasyLight Mitochondriak® devices provide clinical-grade irradiation at 630 nm and 670 nm of visible red light and 810 nm, 830 nm and 850 nm of

invisible near-infrared (NIR) light, or a combination of both. These wavelengths of light provide additive therapeutic effects.

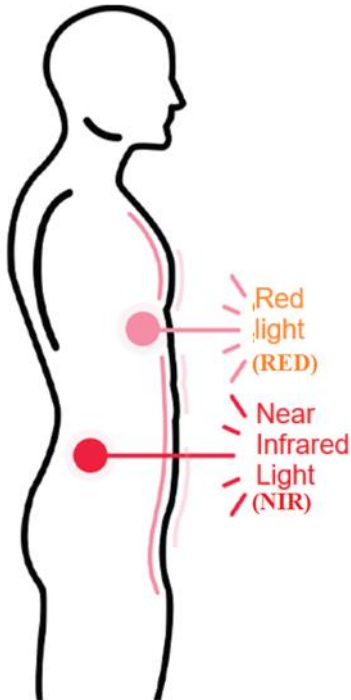Red light 630 and 670 nm:

Most of the red light is absorbed by your skin, leading to increased collagen production and optimal skin health. Red light 670nm also best stimulates ATP production and changes the viscosity of water in mitochondria

Near infrared light 810, 830, 850 nm:

Near infrared light penetrates deeper into your body, improving muscle recovery and reducing joint pain. It also stimulates mitochondria and cellular respiration.

General benefits of therapy:

Red light (RED):

- Supports liver function
- Supports circulation
- Beneficial for bones
- Beneficial for the kidneys
- Skin renewal
- Supports muscle healing
- Helps blood regeneration
- Beneficial for the lungs
- Helps in treating pain
- Beneficial for the heart
- Facilitates wound healing
- Reduces scarring
- Helps improve metabolism
- Enhances the senses

Near infrared light (NIR):

- Pain relief
- Strengthens enzymes
- Relieves allergy symptoms
- Skin renewal
- Improves skin sensitivity to UV
- Helps heal muscles and ligaments
- Helps with sprains
- For faster healing and energy
- Accelerates the healing of burns
- Accelerates skin damage
- Accelerates wound healing
- Strengthens mitochondrial function

Installation and hanging of the device:

5. Screw the bar onto the top of the panel (See pictures below).
6. Screw the suspension rope to the rail with a carabiner.
7. Place the hinge mechanism on the fixed door and use the mechanism to raise the panel to the desired height.
8. The height of the panel can be adjusted (up and down).

Instructions for use:

1. Hang the device in the desired position or place it on solid, level ground.
2. Distance: Apply the light 30-60 centimeters away from the body. The closer you are, the greater the intensity and vice versa.
3. **Duration of use:** Expose each body part to light for 5 to max. 30 minutes, depending on distance. Repeat if necessary.
 - **For treatment of superficial body parts** (energy levels, wrinkles, collagen production, muscle regeneration, pain, testosterone, thyroid, immune function) apply for approx. 5 to 10 minutes (required dose is approx. 1 to 10 Joules)
 - **For treatment of deeper parts of the body** (bone healing, brain performance, fat loss, scarring, sleep quality, muscle regeneration, gut health, hair growth, inflammation) apply for approximately 10 to 30 minutes (required dose is approximately 10 to 60 Joules).
4. Frequency of use: 2 to 14 times per week, as needed or until results are visible
5. Maintaining Results: Consistency is key. Consistency and regularity are especially important to see results.

Frequently asked questions:

Are you just starting out with red light therapy?

It doesn't matter. With every order of the EasyLight Red Light Therapy Panel, you will receive everything you need to access the amazing health benefits of light therapy.

When ordering any device, you will receive a manual, an easylight calculator, a practical pdf and everything else to enjoy your therapy to the fullest. All it takes is to read the entire manual. 😊

How does red light therapy work?

Red and infrared light photons allow our cells to use oxygen efficiently, minimizing the uptake of nitric oxide by cytochrome C oxidase, and thus create sufficient water and ATP. Only red light therapy can penetrate deep into the body and our mitochondria to stimulate ATP production, healing, regeneration, and subsequently help improve your appearance, performance, and overall well-being.

Are you sensitive to light?

Before using the red light therapy panel, we recommend that you check your photosensitivity. To do this, follow these steps:

Expose an area, such as your forearm, 15 cm from the device to the light for 3 minutes. If you experience any redness or discomfort, you may be photosensitized and we recommend that you consult your doctor or healthcare provider before using the device.

What problems can therapy help me with?

Red light therapy is used, for example, to relieve muscle spasms, minor muscle and joint pain, as well as to relieve pain and stiffness associated with arthritis. It is also used during sports or performance to improve regeneration, improve ATP production, or for better overall well-being.

How long should I be enlightened?

The general recommendation for using our panels is: **In general, for optimal results, we recommend using the panel at a distance of 15 to 45 cm from the body and illuminating for approximately 10 to 15 minutes per area** (depending on the size of the panel). **1 to 2 times a day, 5 to 7 days a week.**

You can also easily set the timer on the panel and you don't need to keep track of the time at all. Simply enjoy your therapy, or use it for "meditation" or thinking.

- For treatment of superficial body parts (wrinkles, collagen production, muscle regeneration, pain, testosterone, thyroid, immune function) apply for approx. 5 to 10 minutes (required dose is approx. 1 to 10 Joules)
- For treatment of deeper parts of the body (bones, organs, brain performance, fat loss, scarring, sleep quality, muscle regeneration, gut health and microbiome, hair growth, inflammation) apply for approximately 10 to 30 minutes (required dose is approximately 10 to 60 Joules).

Is consistency really the key?

It's important to remember that, just like with your circadian rhythm and natural exposure to sunlight, it's all about regularity/consistency.

Therefore, we recommend using the red light therapy panel every day or at least 5 times a week. For additional health benefits, including pain relief, aches and tension, even up to 3 treatments per day may be beneficial. Photobiomodulation (PBM) or light therapy is essentially like “food” for your cells and especially the mitochondria.

Therefore, it is important to note that, just like with real food, unfortunately, most people do not get the right amount of healthy light that their body and mitochondria need for optimal well-being. As a result, their body may react with detox-like symptoms. So, if you are experiencing these types of effects, we recommend starting with 1-2 minute treatments per area, slowly working up to 10 minute sessions over 2-3 weeks as your body adjusts to the therapy. As always, if you have any concerns about using photobiomodulation, we recommend consulting your trusted healthcare provider or doctor.

Results and expectations?

The time it takes to see benefits and improvements depends on your health goals. The benefits of light therapy for cellular health have been thoroughly proven in hundreds of clinical studies (many of which can be found on easylight.sk), so trust that your body benefits from each therapy, although the effect may not be immediate like drinking coffee. However, there are many factors that can affect how quickly you feel an improvement in your symptoms. Some people will notice benefits within a few days. But overall improvement in health or skin-related problems may take 8-12 weeks of consistent use. By using a red light therapy panel daily, you can maximize the speed of your body's healing.

Can light hurt eyes?

The typical human eye responds to wavelengths from 400nm to 700nm. The light from this device is visible red or invisible infrared. So your eyes will only see some of the wavelengths from the panel.

Infrared light is beneficial for the eyes as well as other parts of the body. Protective glasses are not required to use with infrared light, but we recommend using them.

The red light is very bright, so it may be uncomfortable to use with your eyes open. Especially if you are sensitive to light or have damaged eyes. Usually, closing your eyes is sufficient to use the light directly on your face. However, if you are unsure or the light bothers you too much, we recommend consulting a specialist or doctor.

What if I overdo it with the enlightenment time?

Since red light therapy does not contain any UV light or produce excessive heat, there are no risks associated with too much red light. However, there is such a thing as "overdosing," and like with anything else (including pure water), if something is good, more is not necessarily better. Overusing red light therapy will most likely simply lead to less positive effects because your cells will be overstimulated.

Therefore, we simply recommend that you follow the instructions in the manual that comes with your device. This way you will know that you are receiving an effective dose.

You can find more answers on our website: www.mitochondriak.com

What's in the package:

1x LED red and infrared panel	1x Manual
1x Protective glasses	1x Hand forging
1x Power cable	1x Door hanger
1x Height-adjustable pulley system	1x Remote control

Maintenance and cleaning instructions:

Maintenance

Our panels do not require any maintenance other than routine cleaning and care.

Cleaning the device

5. Turn off your device by sliding the power switch to the off position on the back of the device.
6. Unplug the power adapter from the electrical outlet.
7. Clean the surface of the device only with a dry cloth or duster and make sure the device is dry before use.
8. Avoid using any harsh cleaners or bleach, or damp cloths, as they may damage your device.

Warnings and Cautions

- Easy Light Mitochondriak® products, intended for light therapy, are not medical devices. The use of our products does not in any way replace a visit to a doctor or consultation with a specialist.
- Do not look directly into the LEDs or their reflection in the mirror as this may cause temporary eye irritation.
- Do not allow children to use the device without supervision.
- This device is an electrical device. To prevent electric shock, observe the following points.
- Do not throw the device into water or liquids as you could receive an electric shock.
- Do not wash electrical parts with water or liquids.
- Do not touch the device directly or touch plugs or switches with wet hands.
- Do not touch the device if your skin is wet or damp.
- Do not pull, carry or lift the appliance by the cable. If the cable is damaged, do not use the appliance as this could result in electric shock.
- Do not use the device if it has been dropped, damaged, or if it has come into contact with water or other liquids.
- Do not use the device on or near heated surfaces.
- Do not use this device in areas where it may be exposed to flammable materials or products, or in places where fumes may be present. Explosion and fire may occur.
- This device can not be used by sleeping or unconscious persons.
- Do not use the device on open wounds as the heat may increase circulation and cause bleeding.
- Do not use this device in combination with ointments and balms that contain heat-generating ingredients, as this may cause skin burns.
- Use of this device in a manner other than that described in this manual will void the warranty and may result in serious injury.

Tips:

- The device is suitable for use indoors and in dry environments only.
- If the panel stops working, disconnect it from the power supply.
- Do not place the panel near a heat source.

Warranty:

We guarantee that this equipment will operate within its specifications for a period of 2 years from the date of delivery.

What does this warranty cover?

This limited warranty covers all defects in material or workmanship under normal use as described above during the warranty period. During the warranty period, our company will repair or replace, free of charge, the product or parts of the products that prove to be defective due to faulty material or workmanship under normal use and maintenance.

What is not covered by the warranty?

The limited warranty does not cover any problem that is caused by damage or the result of improper maintenance or self-modification. It also does not cover damage to the product due to natural disasters such as flooding, falling of the product from a height or other natural disasters, theft or loss of the product.

So if your device stops working, simply write to us.

The device also has a 14-day legal period during which you can return it without giving a reason. However, we give you up to 30 days.

This means that if you write to us within the first 30 days of purchase that you want to return the device, we will refund you in full. Simply send us the device in its original packaging and we will send you your money back!

Try new panel models:
Mitochondriak® 2.0 PULSE
&
Mitochondriak® 3.0
& more

Designed for demanding Mitochondriacs, for athletes, for fitness/wellness centers, but also for the home.

*Let's get closer to nature and sunlight
and discover our full potential!*

Mitochondriak®

Easy**Light** Mitochondriak®


...because Light is shaping Life

www.mitochondriak.com

EasyLight

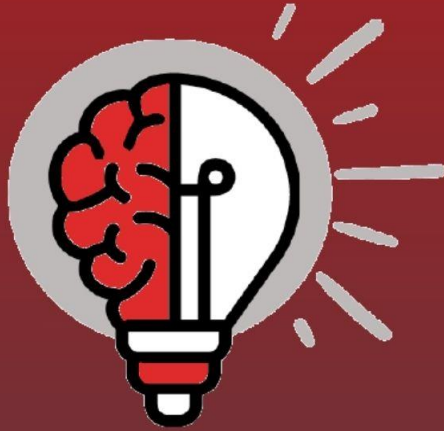Nature is always first
and best answer.

But, what is your second?

Mitochondriak®